Product Information Sheet



Sf-RVN® Insect Cell Line

A rhabdovirus-free cell line improving the safety profile of bioprocesses

Product Description

The Sf-RVN® Insect Cell Line is a rhabdovirus-negative Sf9 cell line that can be grown in both adherent and suspension culture. Sf-RVN® Insect Cell Line has been adapted for growth in chemically defined medium to be used for the production of recombinant proteins, adeno-associated viruses (AAV) and viral-like particles (VLP). The cell line was tested for adventitious agents and banked under GMP conditions. A technical user guide, with detailed protocols for optimal performances, is provided with the cell line. The EX-CELL® CD Insect Cell Medium, recommanded for use, is optimized to get excellent growth and productivity of the Sf-RVN® Insect Cell Line. Combined, these two products form the Sf-RVN® Platform.

Cell Line Origin

Sf-RVN® is an insect cell line originally derived from ovarian tissue of *Spodoptera frugiperda* (*Sf*) *pupae*. The parental line, IPLB-SF-21 (renamed IPLB-SF-21 AE after adaptation into TC-100 medium), was established by J. Vaughn at the Insect Pathology Beltsville in Maryland, USA in the late 1960's. A subclone, designated Sf9, was established in 1983 by G. Smith and C. Cherry. Sf-RVN® was isolated from Sf9 by GlycoBac, LLC.

Legal Information

The Sf-RVN® Insect Cell Line is sold for research use only. A license agreement must be acquired prior to use of the cell line or any products derived from the cell line in a clinical or commercial process. Please contact your local sales representative for more details.

Product Components

Sf-RVN® cells are provided to customers in vials containing 1 mL at 10×10^6 cells/mL. Cells are banked in EX-CELL® CD Insect Cell Medium and 10% DMSO.


Features and Benefits

- Rhabdovirus-negative cell line: Enhanced risk mitigation and biosafety profile
- Adventitious agent tested and GMP banked in EX-CELL® CD Insect Cell Medium
- Full traceability and documentation available for regulatory filings
- EX-CELL® CD Insect Cell Medium: A chemically defined medium optimized to get excellent growth and productivity of the cell line
- Technical user guide with detailed protocols for optimal performances

Sf-RVN® Insect Cell Line Growth

Typical Sf-RVN® culture viabilities are high nineties (95-98%) when cultivated in the EX-CELL® CD Insect Cell Medium with passaging every 3-4 days. A seven day growth assay performed on the cells in duplicate showed day 7 cell densities at 8 x 10^6 cells/mL.

Sf-RVN® Insect Cell Line in EX-CELL® CD Insect Cell Medium Growth Assay

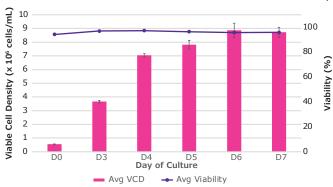



Ordering Information

Product	Cat. No.
Sf-RVN® Insect Cell line	
Sf-RVN® Insect Cell line - 1 Vial	SFRVN-1VL
Related Product	
EX-CELL® CD Insect Cell Medium - Dry Powder	24381C
EX-CELL® CD Insect Cell Medium - Liquid	14380C
Escort™ IV Transfection Reagent	L3287
SyntheChol® NS0 Supplement	S5442

BioReliance® Biosafety Services

Contact **MerckMillipore.com** for additional GMP testing and banking services.

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

To place an order or receive technical assistance in the U.S. and Canada, call toll-free 1-800-645-5476 For other countries across Europe and the world, please visit: **EMDMillipore.com/offices**For Technical Service, please visit: **EMDMillipore.com/techservice**

MilliporeSigma 400 Summit Drive Burlington, MA 01803

EMDMillipore.com



