



GREAT SOLVENTS SMALLER FOOTPRINT

New bio-renewable
Cyrene™ solvent blends
for responsible
scientific research



MilliporeSigma is
the U.S. and Canada
Life Science business
of Merck KGaA,
Darmstadt, Germany.

Sigma-Aldrich®
Lab & Production Materials

A SIMPLE SWITCH. A HUGE DIFFERENCE. **BIORENEWABLE SOLVENTS.**

Dipolar aprotic solvents are regularly used in many research and production applications, from cross-coupling reactions to graphene ink printing. However, petroleum-based solvents, such as DMF (*N,N*-Dimethylformamide) and NMP (1-Methyl-2-pyrrolididinone) are known to have high toxicity.

We make it easy to protect you and the planet with a growing portfolio of safer, greener solvents, like Cyrene™ and 2-Methyltetrahydrofuran (2MeTHF). Now, we have expanded the range to offer you new Cyrene™ blends, and give you even more reasons to switch. Discover how our latest innovations will keep your work performing – sustainably.

Mixed to match your needs

Cyrene™ γ -Valerolactone & Cyrene™ 2-MeTHF blends

Cyrene™ 2-Methyltetrahydrofuran Blend and Cyrene™ γ -Valerolactone Blend combine the benefits of classic Cyrene™ with those of two BioRenewable solvents. While they maintain the excellent efficacy of Cyrene™, the new blends offer considerably lower viscosity. This greatly simplifies their handling and broadens their applicability, leading to even higher product yields in automated processes. Experience greater performance with lower environmental impact by switching from NMP and DMF to our new Cyrene™ blends.

Features & Benefits

- Independently certified as 100% renewable carbon
- Sustainably produced and less hazardous for the environment
- Non-toxic, much safer to handle than petroleum-based aprotic solvents
- Physical properties more analogous to those of NMP and DMF
- Same high reaction yields as traditional solvents
- Optimized for automated processes due to lower viscosity
- Lower boiling point facilitates removal from automated systems





Applications


- HATU amide coupling
- Suzuki–Miyaura coupling
- Sonogashira coupling

Description	Viscosity (cP)	Density (g/mL)	Product No.
Cyrene™ γ-Valerolactone Blend	4.05	1.14	920207
Cyrene™ 2-MeTHF Blend	3.72	1.15	920193
Cyrene™, BioRenewable	10.5	1.25	807796
<i>N,N</i> -Dimethylformamide (DMF)	0.9	0.94	319937
1-Methyl-2-pyrrolidinone (NMP)	1.7	1.03	443778

Read our tech article to learn how the new Cyrene™ blends performed in cross-coupling reactions.
 >> [Link to application note](#)



Suzuki–Miyaura Coupling: Synthesis of 4-Phenyltoluene

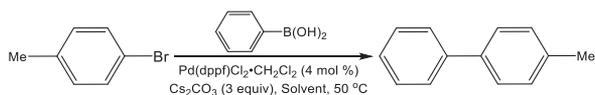

Entry	Solvent	Isolated Yield
1	Cyrene™	85.0%
2	Cyrene™ γ-Valerolactone Blend	99.1%
3	Cyrene™ 2-MeTHF Blend	95.4%
4	<i>N,N</i> -Dimethylformamide (DMF)	99.7%
5	1-Methyl-2-pyrrolidinone (NMP)	97.8%

To discover full range of BioRenewable solvents, visit SigmaAldrich.com/biorenewable



Sigma-Aldrich®

Lab & Production Materials

MilliporeSigma
400 Summit Drive
Burlington, MA 01803

[SigmaAldrich.com](https://www.SigmaAldrich.com)



BIG SCIENCE SMALL FOOTPRINT



[SigmaAldrich.com/
sustainable-chemistry](https://www.SigmaAldrich.com/sustainable-chemistry)

© 2022 Merck KGaA, Darmstadt, Germany and its affiliates. All Rights Reserved.
MilliporeSigma, the vibrant M, Cyrene, and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt,
Germany or its affiliates. All other trademarks are the property of their respective owners.
Detailed information on trademarks is available via publicly accessible resources.

Lit. No. MS_FL8991EN Ver. 1.0
03/2022