


# PrePak L Wall Mounting Bracket

(Cat. No. WBPRPKL01)

## General Description

The PrePak L wall bracket includes the wall bracket, fittings, tubing and pressure gauges needed for the installation and operation of MilliporeSigma's PrePak L pretreatment pack. The PrePak L pack and wall bracket is installed to add additional pretreatment to RiOs™ and Elix® water purification systems. A MilliporeSigma Application Specialist can determine the complete pretreatment required based on feed water contaminant levels.

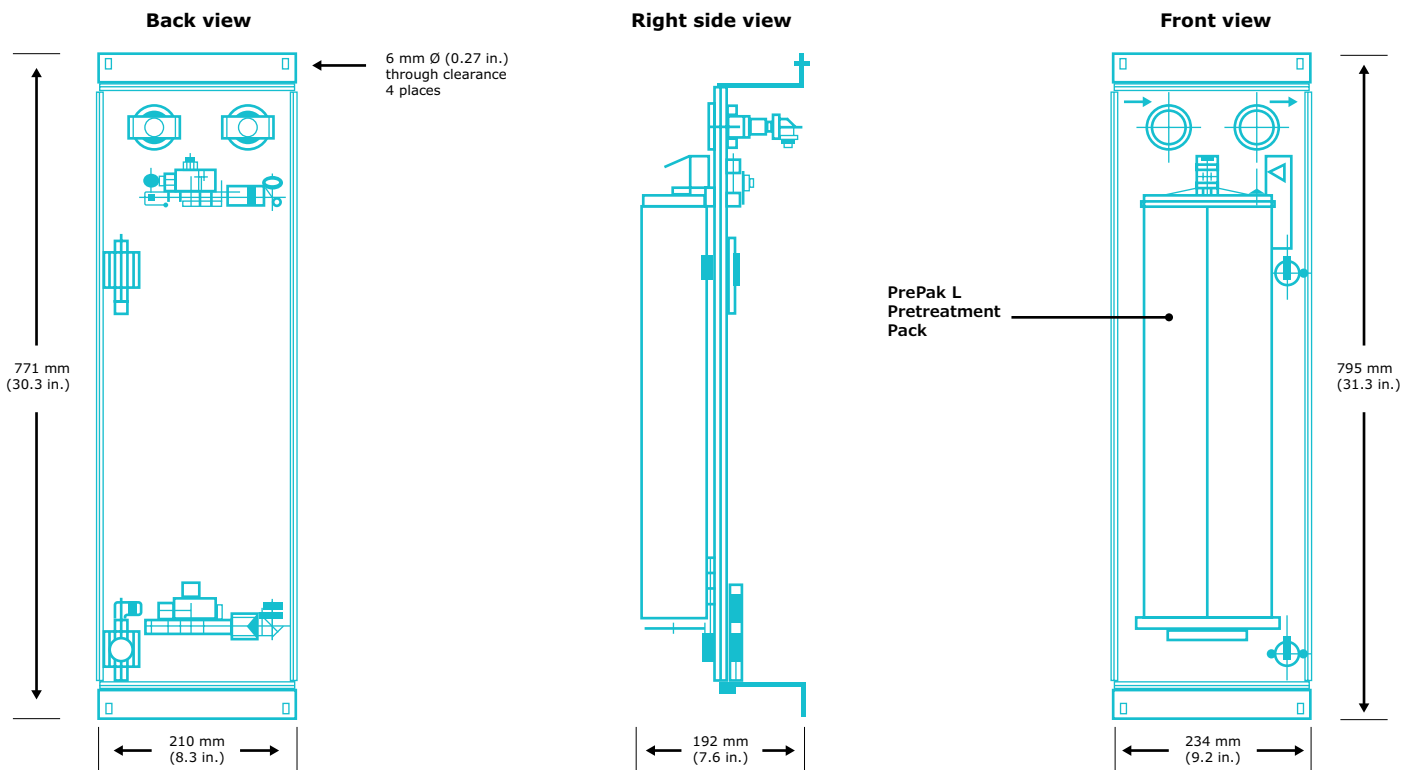

### Hydraulic specifications

Maximum feed pressure 6 bar (86 psi)

Follow RiOs™ and Elix® system specifications for maximum tubing length requirements.

	Connection/component	Description
A	½ in. tube fitting	PrePak L inlet (feed water)
B	½ in. tube fitting	PrePak L outlet (to water system)
C	Panel mount pressure gauge	PrePak L pressure gauges (inlet and outlet)
D	Panel mount ball valve	PrePak L purge-drain ball valves
E	8 mm tube fitting	PrePak L purge-drain valve
F	Locking mechanism	PrePak L lock/un-lock mechanism

### Installation space requirements

- Additional space will be required for inlet and outlet tubing and purge drain connections.
- The PrePak L wall bracket can be wall-mounted near the inlet of MilliporeSigma RiOs™ and Elix® systems. Use the dimensions from the Back View diagram above for mounting hole locations.
- The PrePak L wall bracket can also be easily mounted to the system cart (MilliporeSigma catalog number ZSCART001). Contact MilliporeSigma for specifications.

## For More Information

Please visit our website  
[EMDMillipore.com/labwater](http://EMDMillipore.com/labwater)

MilliporeSigma  
290 Concord Road  
Billerica, MA 01821

